



仪器方法:

液相色谱仪

色谱柱: 两通

流动相:

Time(min)	A: 水	B: 乙腈
0	20%	80%
2	20%	80%

流速: 0.6mL/min

柱温: 30°C

检测波长: DAD320nm

溶剂: 甲醇+水

MS

Ion source: ESI

MS: Scan

Mass(m/z): 50-400


Cone Voltage(V): 20

Ionization Mode: ES+/ES-

Time (min) : 0min-2min



上海源叶生物科技有限公司
Shanghai yuanye Bio-Technology Co., Ltd
电话: 400-666-5481 传真: 021-55068248
网址: www.shyuanye.com
邮箱: 3008007435@qq.com

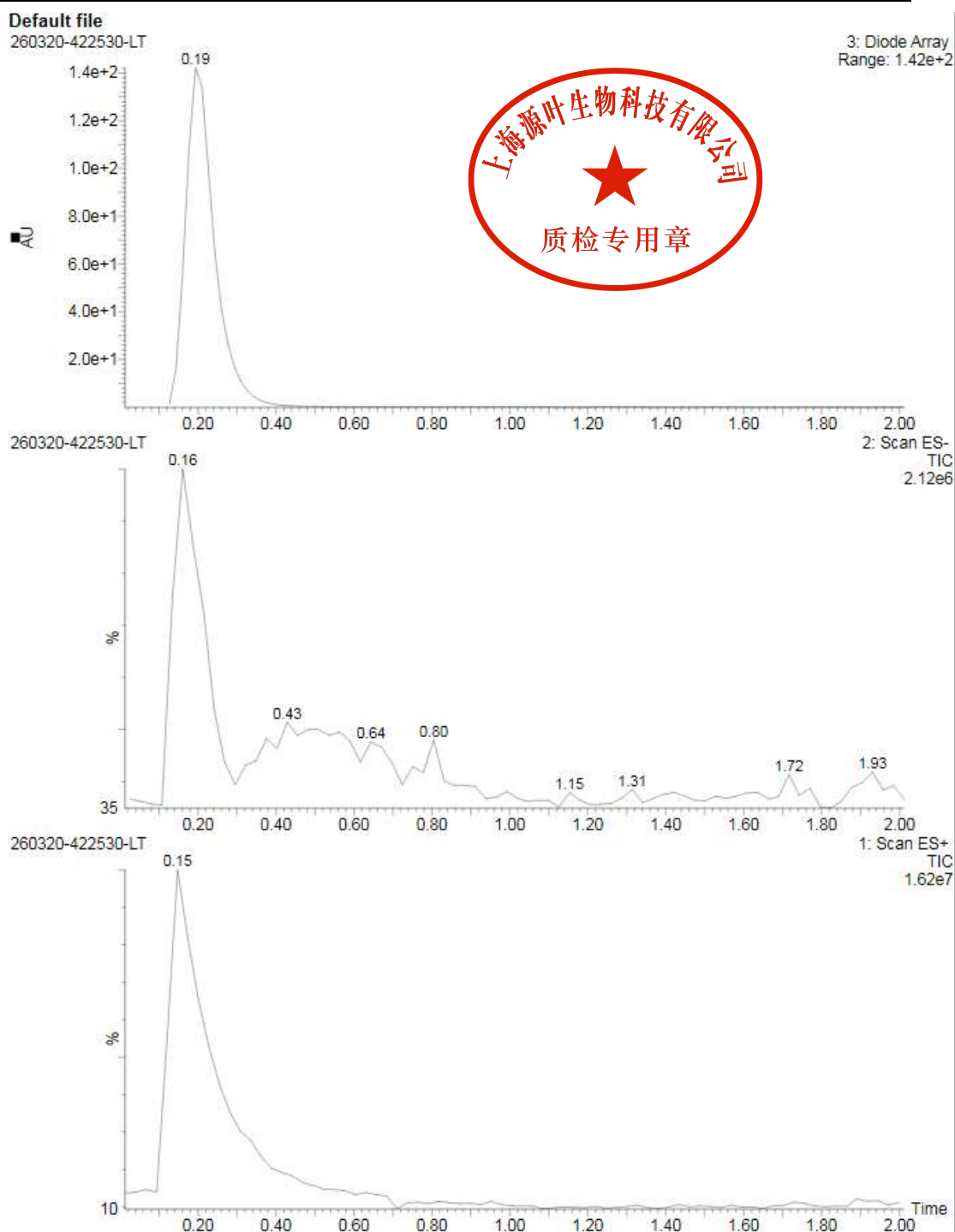



上海源叶生物科技有限公司
Shanghai yuanye Bio-Technology Co., Ltd
电话: 400-666-5481 传真: 021-55068248
网址: www.shyuanye.com
邮箱: 3008007435@qq.com

